
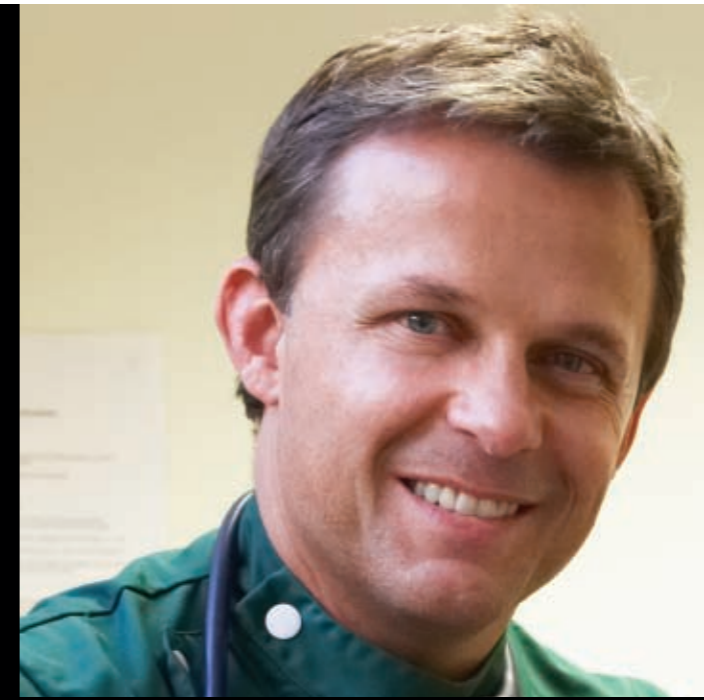


Established in the market  
place since 1989

Centurion Scientific Ltd

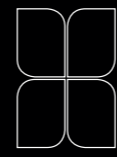

Distributed by


Manufactured by:

Centurion Scientific Ltd

email: [info@centurionscientific.co.uk](mailto:info@centurionscientific.co.uk)  
[www.centurionscientific.co.uk](http://www.centurionscientific.co.uk)



Pro-Vet

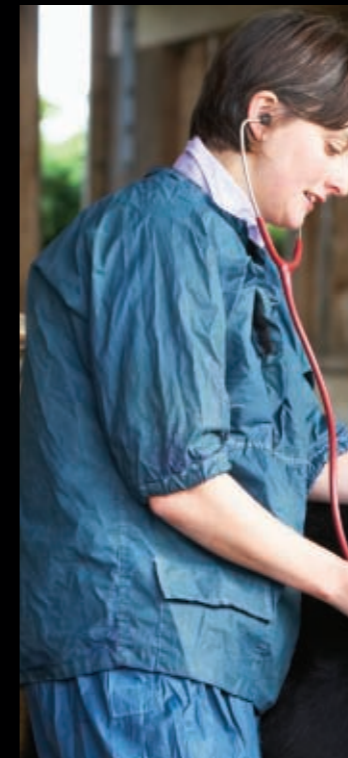

Designed by  
Centrifuge Engineers  
and Veterinary  
practitioners...

...together combining their experience to bring you the ultimate Veterinary Centrifuge's

The Pro-Vet HE offers small animal use from Capillary tubes (75mm) to 2ml

The Pro-Vet Multi offers the same as the HE, but adding a swing out rotor to give increased sample capacity for larger animals, up to 100ml.

## Pro-Vet HE

The dual purpose veterinary centrifuge that separate's blood, urine, haematocrit and many other routine samples.

The standard, inclusive Combination rotor is simple simplicity not only allowing capillary tube analysis (PCV) but offering extra capacity of 2ml tubes , for larger animal analysis. Thus saving you time and money by not requiring further rotors and saving time in changing them.



## Pro-Vet HE

With Bright blue LED display

<b>Speed</b>	<b>500-12000Rpm</b> in 10rpm steps
<b>Max Rcf</b>	<b>13500G</b>
<b>Timer</b>	<b>0-99mins</b> and constant. Pulse 0-90 secs.
<b>Dims HWD</b>	<b>23x23x35cms</b>
<b>Weight</b>	<b>12 Kilos</b> without rotor

### Features to both models

- LED Display
- 10 acceleration rates
- 10 deceleration rates
- 10 programme memory
- Run in speed or Rcf (G force)

## Pro-Vet Multi

With Bright blue LED display

<b>Speed</b>	<b>500-12000Rpm</b> in 10rpm steps
<b>Max Rcf</b>	<b>13500G</b>
<b>Timer</b>	<b>0-99mins</b> and constant. Pulse 0-90 secs.
<b>Dims HWD</b>	<b>27x33x46cms</b>
<b>Weight</b>	<b>16.5 Kilos</b> without rotor

### Pro-Vet Multi only

- Rotor recognition system

### Standard features

- Stainless steel bowl
- Alloy & steel main frame with composite panels
- World brand leading industrial grade inverter
- High quality European brushless motor
- Full service access via coded system

### Safety features to both

- Multi Point Lid locking system
- Imbalance vibration detector with auto reset
- Lid lock detector
- Overspeed sensor
- Set Inverter values
- Barrier Ring
- Motor overheat sensor


### Advantages

- Cleanliness & rust free
- Strong, light & quiet
- Reliability & accuracy
- Quiet, cool & reliable
- Easy servicing

### Advantages


- Best safety for lid
- Easy "eye" volume loading
- Unable to run with an open lid
- Safe detection of speed
- Electronic safety of speed
- Metal protection of chamber
- Motor & samples protection

### Complete with the following

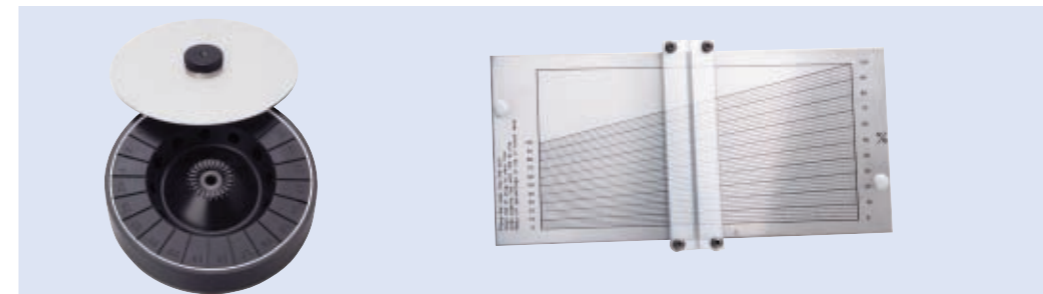

BR5401 Combination rotor

8003 PCV Reader

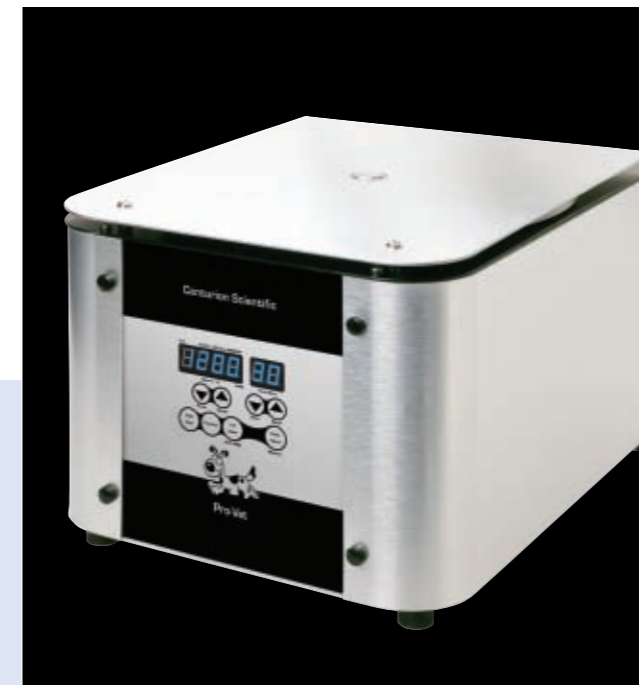

## Pro-Vet Multi

Equine, Zoological, Larger animals. This machine offers fantastic flexibility by using the combination rotor and a swing out rotor for those larger samples from up to 15ml of blood to 50-100ml of sperm.

All in complete safety as each rotor has a recognition system matched to the centrifuge. Just change the rotor that's all you have to do.

Having one machine that can do both large and small animal's means huge space and cost savings.

### Complete with the following



BRK5401 Combination rotor

8003 PCV Reader

BRK5505 Rotor

Adaptors (1set)